


preVent™ Filter

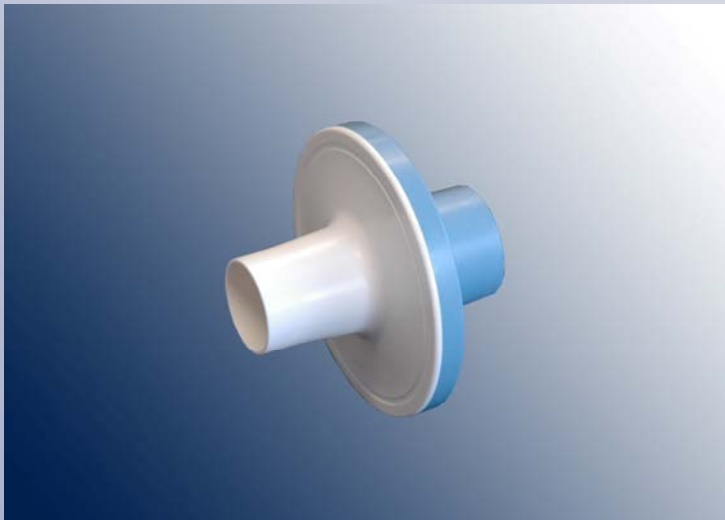

All specifications subject to change without notice. Products may vary from that illustrated. Please contact your Medical Graphics representative for latest information, pricing and product availability.

©Copyright, Medical Graphics Corporation, 2006. MedGraphics is a registered trademark of Medical Graphics Corporation.

Part #060027-001 REV B

Introducing the MedGraphics preVent Filter. The custom filter for your MedGraphics product.

The MedGraphics preVent Filter meets or exceeds the industry standards for filtering in pulmonary function testing. The MedGraphics preVent Filter is specially designed to fit on the MedGraphics preVent pneumotach. The preVent Filter negates the need to use a coupler which is needed with other standard filters. The MedGraphics preVent Filter assures a secure fit with no leaks. It is the elegant solution for your infection control needs with your MedGraphics system. The preVent Filter uses a medium with high quality electrostatic properties that traps bacteria and viruses, preventing contamination to both the patient and equipment.

Features

- No need for coupling devices for your MedGraphics System
- Leaks are no longer a concern with the preVent Filter
- Low resistance
- Minimal deadspace
- Meets or exceeds ATS guidelines

MedGraphics preVent Filter Specifications

Bacterial Efficiency

- 99.99%

Resistance

- 0.23 cm H₂O/L/sec approx.

Viral Filtration Efficiency

- 99.99%

Materials

- K resin case
- Electrostat filtrate medium

Dead Space

- Approximately 35cc

MedGraphics Order Information

Part Numbers

- Filter 536713-001
- Kraton Mouthpiece 536101-001

Order Contact Information

- Fax: 651-379-8333
- Email: filters@medgraphics.com

Call your MedGraphics supply representative today at 1-800-655-1133 to discuss special prices on quantity purchases.